
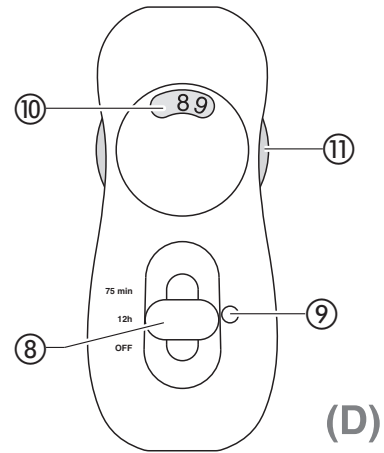
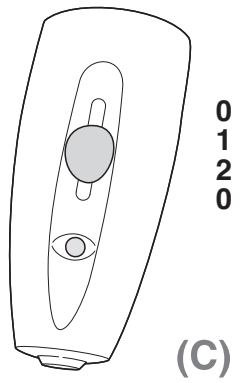
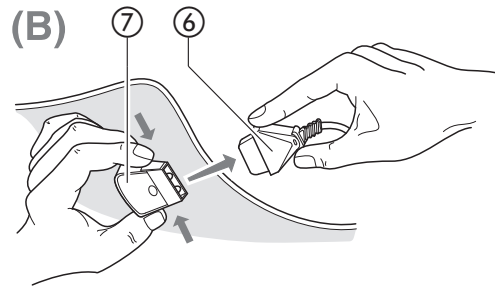
*Istruzioni ed Avvertenze
Instructions and Warnings
Instructions et avertissements
Anleitungen und Hinweise
Instrucciones y Advertencias
Instruções e Advertências
Navodila in opozorila
Οδηγίες και Προειδοποιήσεις
Bilgiler ve Uyarılar
Návod a upozornění
Utasítások és figyelmeztetések
Upute i napomene
Инструкции и рекомендации*

SCALDALETTO

HEATED UNDERBLANKET
CHAUFFE-LIT
BETTWÄRMER
CALIENTACAMAS
CAPA AQUECEDORA
POSTELJNI GRELEC
ΘΕΡΜΑΙΝΟΜΕΝΟ ΚΑΛΥΜΜΑ
ΣΤΡΩΜΑΤΟΣ
YATAK ISITICISI
VYHŘÍVACÍ PROSTĚRADLO
FÚTÓTAKARÓ
GRIJAČ POSTELJE
ГРЕЛКА "НАМАТРАСНИК"

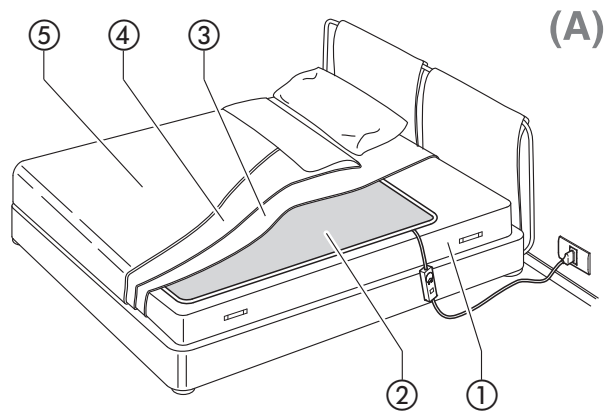ATTENZIONE

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WARNING

INSTRUCTIONS AND WARNINGS FOR CORRECT AND SAFE USE

IMPORTANT INSTRUCTIONS. KEEP FOR FUTURE REFERENCE

THEY SUPPLY IMPORTANT INSTRUCTIONS REGARDING OPERATION AND SAFETY DURING PRODUCT INSTALLATION, USE AND MAINTENANCE.

GENERAL WARNINGS

1. After unpacking, make sure that the appliance is intact with no visible signs of damage caused during transport. Should you have any doubts, do not use the appliance and take it to an authorised IMETEC Service Centre.
2. Before connecting the appliance, check that the rating plate specifications correspond to the electrical mains supply. The rating is located on the label stitched onto the appliance.
3. Do not use adapters, multiple sockets and/or extension cords. If the appliance's plug and the socket are incompatible, have the socket replaced with a suitable type by qualified staff.
4. This appliance must **only** be used as an electric underblanket as indicated in the instructions. Any other use is to be considered hazardous. The manufacturer cannot be held responsible for any damage due to improper or incorrect use.
5. The use of any electric appliance implies the observance of some fundamental rules. Specifically:
 - never touch the appliance with hands or feet that are wet or damp
 - do not tug the power cord or the appliance to pull the plug out of the socket
 - do not expose the appliance to the elements (rain, sun, etc.)
6. Before cleaning or servicing the appliance, disconnect it from the mains by pulling out the plug.
7. Should the electric underblanket or protection fuse break down and/or fail to work properly, turn off the appliance and do not tamper with it. If it needs repairing, take it to an authorised IMETEC Service Centre.
8. To avoid dangerous overheating, we recommend unwinding the power cord completely during use and disconnecting the plug from the socket when the appliance is not in use.
9. If the control power cord is damaged, the control must be replaced by the manufacturer or an authorised IMETEC Service Centre to prevent any risks.
10. Should you decide not to use the appliance any more, we suggest making it inoperative by cutting the power cord after removing the plug from the socket. We also recommend making all potentially dangerous parts of the appliance harmless, particularly for children who might be tempted to play with it.

11. Use the electric underblanket only between the mattress and the bottom sheet, never over persons in bed.
12. Do not use the heated underblanket with invalids, infants or people insensitive to heat. The heated underblanket must not be used by children if the control devices have not been preset by the parents or guardian or if the child has not been adequately taught how to use these devices safely. Children must be supervised to make sure they do not play with the appliance. Wearers of pacemakers are advised not to use the heated underblanket for the whole night. The heated underblanket may be used to preheat the bed, but must be switched off or disconnected from the electricity power supply before lying down.
13. Do not use on an adjustable bed: if you do, make sure that the blanket and the power cord cannot get caught or stuck.
14. Do not tuck in.
15. Avoid creasing the appliance.
16. This product is not designed for medical use in hospitals.
17. The appliance must only be used with the control types indicated on the appliance's label.
18. Do not position the control underneath the mattress or pillow.
19. Frequently check the appliance to see if there are any signs of wear and deterioration. If you find any, or if the appliance has been used incorrectly, return it to the supplier before using it again.
20. The fabric could undergo natural yellowing over the course of time. This occurrence does not affect the safety and/or performance of the product itself in any way.
21. This appliance should not be used by children or persons with partially reduced physical, sensorial or mental faculties, and in any case not without specific knowledge and experience, or without the supervision or instructions for use provided by a responsible person.
22. Children must not be allowed to play with the appliance.
23. When storing the appliance, allow it to cool down before folding

This product complies with EU Directive 2002/96/EC.

The crossed bin symbol on the appliance indicates that the product, at the end of its life span, must be disposed of separately from domestic waste, either by taking it to a separate waste disposal site for electric and electronic appliances or by returning it to your dealer when you buy another similar appliance. The user is responsible for taking the appliance to a special waste disposal site at the end of its life span. If the disused appliance is collected correctly as separate waste, it can be recycled, treated and disposed of ecologically; this avoids a negative impact on both the environment and the health and contributes towards the recycling of the product's materials. For further information regarding the waste disposal services available, contact your local waste disposal agency or the shop where you bought the appliance.

**DO NOT USE IF FOLDED
OR CRUMPLED**

**DO NOT INSERT
NEEDLES OR PINS**

IMETEC thanks you for choosing this product. We are certain that you will appreciate the quality, reliability and regard for the consumer with which it was designed and produced.

HOW TO POSITION THE ELECTRIC UNDERBLANKET (Fig. A) (All models)

Spread the electric underblanket ② directly above the mattress ① (do not position it underneath the pillow) with the side containing the plug ⑦ turned towards the mattress. Make sure that it is smoothly spread out. Lay the bottom sheet ③ on top of the electric underblanket, then finish making the bed as normal with another sheet ④ and bedcovers ⑤.

HOW TO CONNECT OR REMOVE THE TEMPERATURE CONTROL DEVICE (Fig. B) (All models)

THE ELECTRIC UNDERBLANKET CAN ONLY BE USED WITH THE TYPE OF CONTROL OR SWITCH SUPPLIED, AS SPECIFIED ON THE LABEL STITCHED ONTO THE APPLIANCE.

To connect the temperature control to the electric underblanket insert the control device socket ⑥ in the plug ⑦.

For the double blanket, repeat this operation for the other socket. For easier bed making, we suggest you remove the socket ⑥ from the electric underblanket by pressing the clips on either side of the plug ⑦.

HOW TO SWITCH ON THE ELECTRIC UNDERBLANKET AND CHOOSE THE RIGHT TEMPERATURE

Important: Do not place the control device underneath the pillow or bedcovers when it is turned on. *Before reading any further, check the type of control supplied with your electric underblanket.*

Dual temperature control (Fig.C)

Plug into the mains socket. Switch the electric underblanket on at least one hour before going to bed setting the temperature control switch to level 2 (maximum temperature). After about one hour, the bed will be pleasantly warm.

If the room is very cold, turn on the electric underblanket two or three hours before going to bed.

If you wish to keep the electric underblanket warm all night, turn the control to setting 1 (minimum temperature) when you go to bed. Do not use setting 2 (maximum temperature). If you want to turn the electric underblanket off, turn the control to setting 0 and remove the plug from the socket.

Control with automatic switch-off timer (Fig. D)

The control with a timer has two automatic switching off phases: after 75 minutes or after 12 hours. To switch the electric underblanket on, put the plug into the mains socket. Set the switch ⑧ to the required position according to whether automatic switch-off is to take place after 75 minutes or 12 hours. Activation of the electric underblanket is shown by the orange LED ⑩ lighting up. Once the required time has been selected the electric underblanket will operate until the time elapses, after which it will switch off. Switching off is signalled by the indicator light ⑨ beside the control key blinking constantly. To set the required temperature (Fig. D), turn the temperature regulation knob ⑪.

During the bed pre-heating phase it is recommended to use the maximum temperatures (8/9) while you are recommended to lower the temperature to between 1 and 4 when going to bed, as shown in table A.

Table A

TEMPERATURE	
1-4	Night use
5-7	Pre-heating - medium temperatures
8-9	Pre-heating - high temperatures

The bed will be pleasantly warm after about 10 minutes, but 20/30 minutes are needed before the maximum temperature selected is reached. If you want to turn the electric underblanket off, position the key ⑧ on OFF and **remove the plug from the socket.** If you want to reprogram

automatic switch-off again set the key ⑧ to OFF, wait for a few seconds and set it to the required position again. You are recommended to switch the electric underblanket off the next day when you get out of bed.

INSTRUCTIONS AND WARNINGS FOR WASHING AND DRYING

DO NOT DRY
CLEAN

HAND WASH
(40°)


DO NOT USE
CHLORINE

DO NOT
WRING

DO NOT
TUMBLE DRY

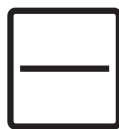DRY FLAT

DO NOT IRON

The electric underblanket can be washed, since all the electric parts and the plugs have watertight seals. We recommend hand washing all wool models in lukewarm water (40°). Only the polyester models can be machine washed using the delicate program (40°). The electric underblanket should only be washed if it is particularly dirty; too many washes can wear out the product (in any case, the electric underblanket should not be washed more than two/three times a year).

Hand wash as follows:

- DETACH THE CONTROL DEVICE FROM THE ELECTRIC UNDERBLANKET - (Fig. B)
- Use a neutral detergent dissolved in lukewarm water (max. 40°C).
- Soak the electric underblanket for 15 minutes.
- Rinse thoroughly in cool water without wringing or squeezing: just shake it to remove the water.
- Lay the electric underblanket out flat and leave it to dry naturally outside for one or two days until it is completely dry.
- DO NOT USE THE ELECTRIC UNDERBLANKET IF IT IS DAMP.

HOW TO STORE THE ELECTRIC UNDERBLANKET

The best way to store your electric underblanket is to leave it on the bed. In any case, store it in a dry place and do not rest heavy objects on it which could cause heavy creasing.

REGULAR CHECKS

Check the appliance and power cord frequently for any signs of wear or deterioration. Should there be any such signs, or if the product fails to work properly, take it to an authorised IMETEC Service Centre.

In any case you are recommended to have your electric underblanket checked by an authorised IMETEC Service Centre at least once every 3 years.

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