

Instruzioni per l'uso
Operating instructions
Instrucciones para el uso
Инструкции по применению
Uputa za uporabu
Návod za uporabu
Οδηγίες χρήσης
Instruções para o uso
Návod k pouziti

BILANCIA PESAPERSONE

BILANCIA PESAPERSONE
BATHROOM SCALE
BALANZA PESA PERSONAS
НАПОЛЬНЫЕ ВЕСЫ
OSOBNA TEHTINA
OSOBNA VAGA
ΖΥΓΑΡΙΑ ΠΑΤΗ ΖΥΓΙΣΤΗΡΟΣ
BALANCA PARA PESAR PESSOAS
OSOBNÍ VAHA

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Type M4301

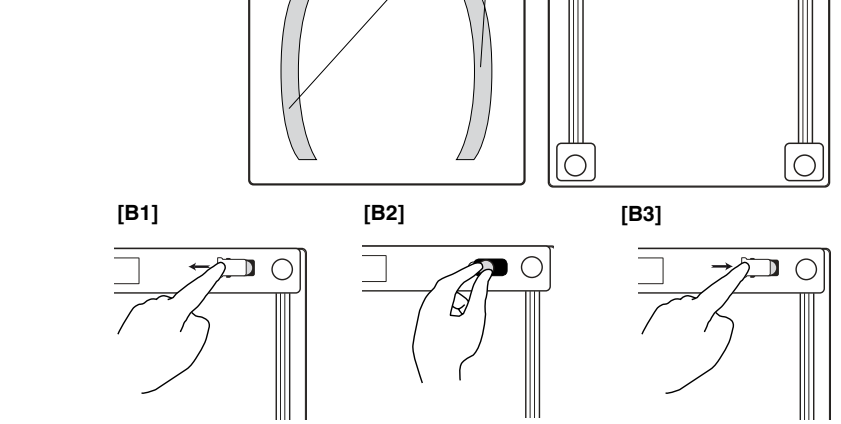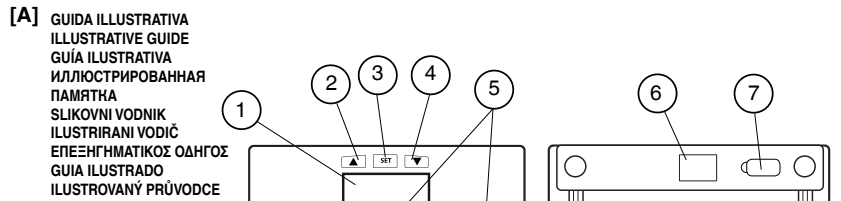[C] LEGENDA DISPLAY LCD (1) / LCD DISPLAY KEY (1) / LEYENDA DE LA PANTALLA LCD (1) / LEGENDA УСЛОВНЫХ ОБОЗНАЧЕНИЙ НА ЖИ ДИСПЛЕЕ (1) / LEGENDA LCD PRIZKAZALNIKA (1) / TUMAS ZASLONA LCD (1) / ΛΕΓΑΝΤΑ ΟΘΟΝΗΣ LCD (1) / LEGENDA ECRA LCD (1) / VYSVĚTLIVKY LCD DISPLEJE (1)

Age	Fat	Water	Muscle	cm	st	lb
18.00	%	%	%	kg	kg	kg

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ENGLISH

OPERATING INSTRUCTIONS FOR THE BATHROOM SCALE

Dear customer, IMETEC thanks you for having purchased this product. We are certain you will appreciate the quality and reliability of this appliance, designed and manufactured with customer satisfaction in mind.

ATTENTION! Instructions and warnings for safe use

Before using this appliance, carefully read the instructions and in particular the safety warnings, which must be complied with. Keep this manual and illustrative guide with the appliance for future consultation. Should you pass the appliance to another user, make sure to also include this documentation.

SAFETY WARNINGS

After having removed the appliance from its packaging, verify its integrity both in the drawing and any presence of damage caused by transport. If in doubt, do not use the appliance and contact an authorised service centre.

This appliance must be used only for the purpose for which it was designed, namely as bathroom scale for home use. Any other use is considered non-compliant and therefore dangerous.

The appliance can be used by children aged from 8 years and above and provided that they are supervised by an adult.

The appliance is faulty or malfunctions, switch it off and do not tamper with it. Any repairs must only be carried out by an authorised technical service centre.

The bathroom scale is not to be used for public or medical use or commercial transactions.

Protect the appliance from cosmetics, humidity, dust, chemicals, strong temperature changes and heat sources (e.g. stoves, heaters).

Place the scale on flat and stable surfaces only. Do not step on the bathroom scale with wet feet if the surface of the scale is wet slipping hazard!

Never leave your child unattended on the scale!

IMPORTANT! The bathroom scale must not be used by persons with a pacemaker or other medical devices.

If the bathroom scale is not used for a long period of time, remove the batteries from the battery compartment.

Replace a battery only with one of the same type.

Make sure to insert the batteries correctly, referring to the symbols for positioning the positive (+) and negative (-) of each battery.

Store the batteries at room temperature in a cool and dry place. Remove the batteries from devices that are not used for long periods of time.

DO NOT charge non-rechargeable batteries.

DO NOT change the non-calibrable batteries with equipment that is not specifically indicated for this purpose or in a different way from that described in the instructions.

NEVER expose the batteries to excessive heat, such as that caused by sunlight, fire or similar sources which pose a high risk of acid leaking from the batteries! Failure to comply with these instructions can damage the batteries and cause them to explode.

DO NOT disassemble or short circuit the batteries or throw them in fire.

DO NOT expose the batteries to excessive heat, such as that caused by sunlight, fire or similar sources which pose a high risk of acid leaking from the batteries! Failure to comply with these instructions can damage the batteries and cause them to explode.

DO NOT disassemble or short circuit the batteries or throw them in fire.

DO NOT expose the batteries to excessive heat, such as that caused by sunlight, fire or similar sources which pose a high risk of acid leaking from the batteries! Failure to comply with these instructions can damage the batteries and cause them to explode.

The acid contained inside the batteries is corrosive. Prevent any contact with skin, eyes or clothes.

Remove the batteries from the appliance before disposing of it.

Dispose of the batteries at a specific collection point.

SYMBOLS

DESCRIPTION OF THE APPLIANCE AND ITS ACCESSORIES

Refer to the illustrative Guide [A] to check the content of the package.

1. LCD display
2. Function key ▲
3. SET function key ▼
4. Function key ▼

USEFUL TIPS

Always weigh yourself in the same conditions, preferably after waking up, in order not to consider weight fluctuations that take place during the day related to food, digestion or physical activity of each individual. Changes in weight may be detected from one day to another.

To use the functions to measure your body characteristics, stand on the bathroom scale with bare feet, taking care to position yourself correctly on the electrodes (5).

During this operation remain perfectly still on the scale.

The display time depends on your ability to stabilise on the bathroom scale and the weight itself; generally, the heavier the weight, the longer the display time.

An incorrect measurement may be displayed in the following cases:

• Persons with a temperature.

• At the end of the age of 10 (the appliance is not calibrated to operate below this age).

• Pregnant women (however, there is no danger).

• Persons undergoing dialysis, persons with limb oedemas and persons suffering from diabetes.

When reviewing one's physical performance, women should consider that during their menstruation they may experience water retention, in which case this may affect the body weight by up to 1 kilogram.

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